ReplengenTM

Hyaluronic Acid Dermal Filler

Contents

1. Hyaluronic acid Filler ?
2. Advantages of Only Replengen™
3. Replengen™ Filler
4. Replengen™ Smart Portable Injector

Hyaluronic acid Filler · HA 필러

- The filler is a syringe type medical device. When injected into the body, it temporarily replenishes the skin tissue to improve facial and body wrinkles and restore the volume.

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- The filler is deep wrinkles, low foreheads and noses and others fills lower depressed areas that can not be solved with Botox, and the procedure time is short within 15 minutes and the effect can be checked immediately after the procedure

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Hyaluronic acid · 히알루론산

- Hyaluronic acid is an in vivo substance widely distributed in human body, chicken crest, skin tissue, placenta, shark cartilage and is harmless to human body.

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
- Hyaluronic acid has excellent biocompatibility and very high—viscoelasticity even when it is transplanted or injected into living body is used for various purposes because it has moisture retention, antibacterial, anti—inflammation and lubrication effect.

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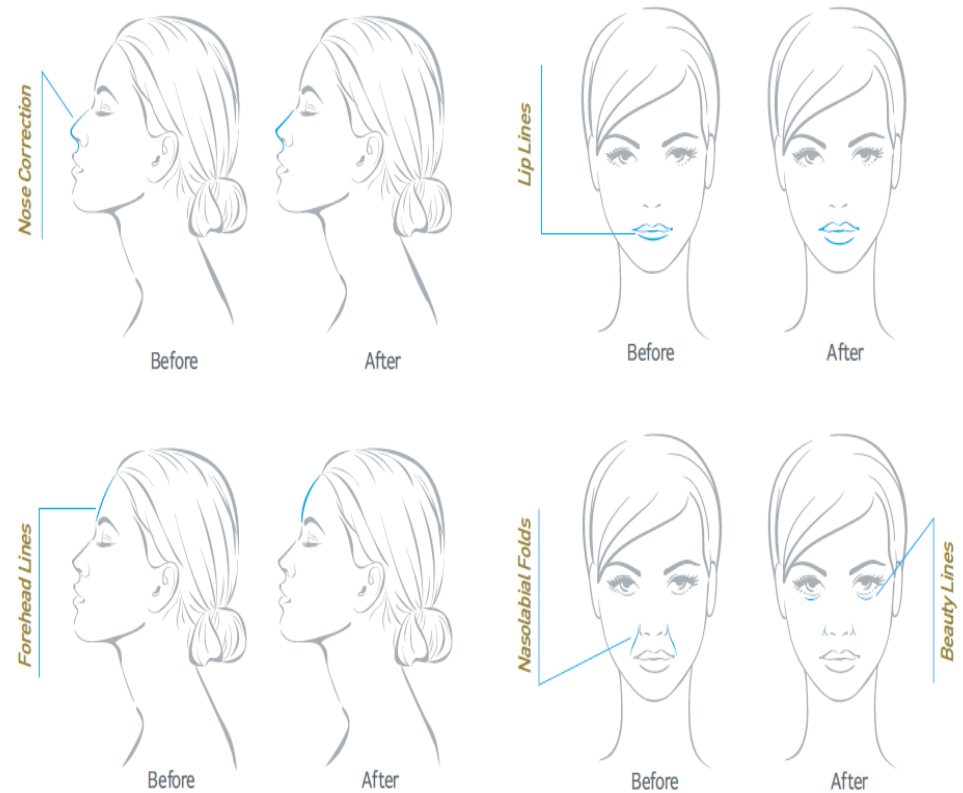2. Advantages of Only Replengen™

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- Glenzen's GI (Grinding Irregular) technology is used to fine-tune the gel particles of the filler to achieve precise pressure output, ensuring sooth injection and stable operation.

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- Replengen™ used the highest quality hyaluronic acid(Made in Japan / Shiseido company, Ltd) and safety was secured by completely eliminating the risk factors through real-time monitoring of harmful substances such as BDDE, endotoxin, heavy metals and microbial contamination that can occur in the manufacturing process.

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Excellent Quality · 우수한 품질

- The Replengen™ filler, which was created by research and development of the best R&D manpower in the industry, introduction the latest experimental equipment and has been proved through strict quality control and guarantee system and maximized the effect of wrinkle improvement and volume recovery after the procedure.

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Line-up 1.1mL

Replegen™ **Fine** With Lidocaine

Used to treat shallow wrinkles such as eye wrinkles, neck wrinkles and fine wrinkles.

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Replegen™ **Deep** With Lidocaine

Used to treat moderate or deep wrinkles such as nasolabial fold, or enlargement of the face, chin, lips.

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Replegen™ **Volume** With Lidocaine

Excellent for creating and maintaining structures with high viscoelasticity and used to improve deep wrinkles and facial contours.

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Replegen™ **Implant S** With Lidocaine

The most viscous product used to improve deep wrinkles and facial contour.

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Line-up 10.0mL

Replegen™ 20

Used to improve the body's shallow wrinkles and body contouring.

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Replegen™ 40

Used to improve the body's deep wrinkles and body contouring.

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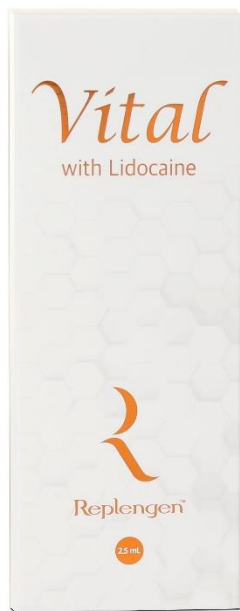
Specification

Product Name	Replengen™ Vital With Lidocaine
Gel Type	Monophasic
Crosslinked HA	20mg/mL
Lidocaine HCl	3mg/mL
Viscoelasticity	10 ~ 30 Pa
Volume	2.5mL
Packing	2 syringe/ carton
Duration	2~4months
Crosslinking Agent	BDDE (1,4-ButaneDiol Diglycidyl Ether)
Injection area	Shallow wrinkles
Injection Depth	Epidermis
Storage	2~25℃
Expiry date	24 months from the date of manufacture

3. Replegen™ Filler

Safety Data

Test Item		Criteria	Result
Appearance		Colorless transparent, Viscous liquid with no foreign material	Conformity
pH		6.8~7.5	7.01
Actual capacity		More than 1.1mL	1.1mL
Osmolality		200~400 mOsm/kg	302 mOsm/kg
Assay	1) Sodium Hyaluronate	Sodium Hyaluronate 90.0 ~ 110.0%	101.1%
	2) Lidocaine HCl	Lidocaine HCl 0.27 ~ 0.33% (2.7 ~ 3.3mg/ML)	0.31%
Endotoxin		< 12.5 EU/ML	< 0.05 EU/mL
BDDE residue		< 2.0 ppm	0 ppm
Heavy metal		< 10.0 ppm	0 ppm
Sterility test		sterile	Conformity
Packaging		Suitable for criteria	Conformity

4. Replengen™ Smart Portable Injector


Product Information

<u>Product name</u>	▪ Replengen Smart Portable Injector
<u>Purpose of use</u>	▪ Injection of drug · 약물주입기
<u>Manufacturer</u>	▪ GST Co., Ltd. · (주)그랑에스엔티
<u>Manufacture and sales</u>	▪ Reanzen Co., Ltd. · (주)리엔젠

4. Replengen™ Smart Portable Injector

9-pin exclusive multi needle · 전용 9Pin Needle

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
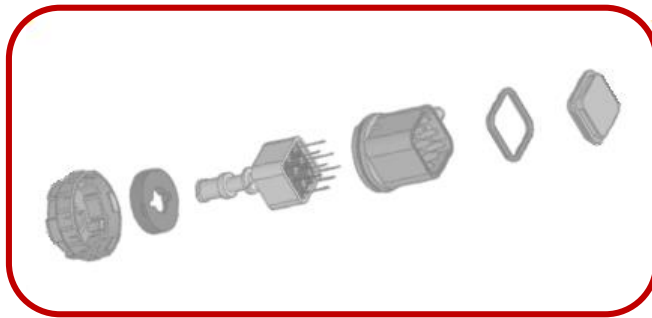
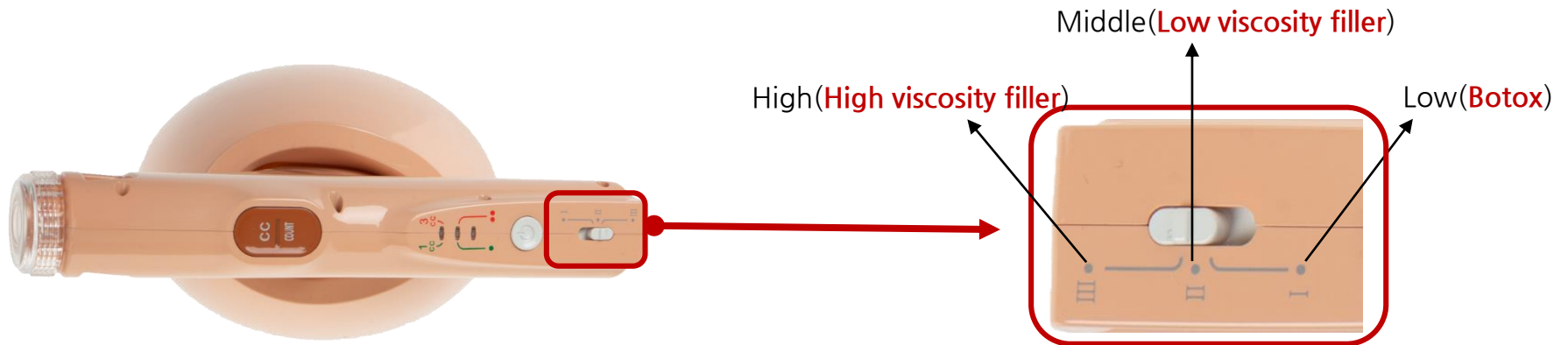
Suction function · 흡입 기능

- It works when the shot switch is pressed.
샷 버튼 누를때 작동
- It can be evenly injected throughout the face.
얼굴 전체에 고르게 주입

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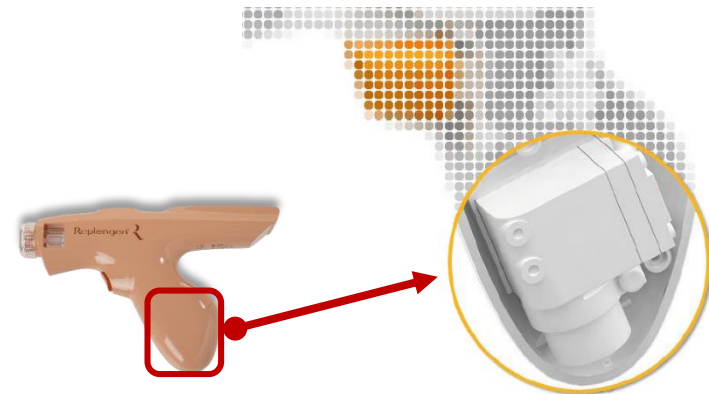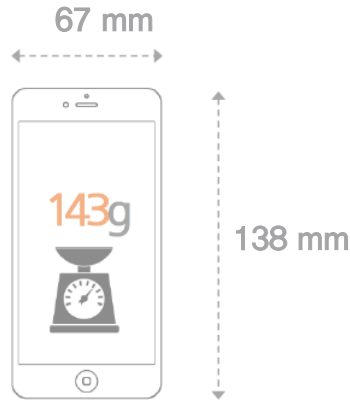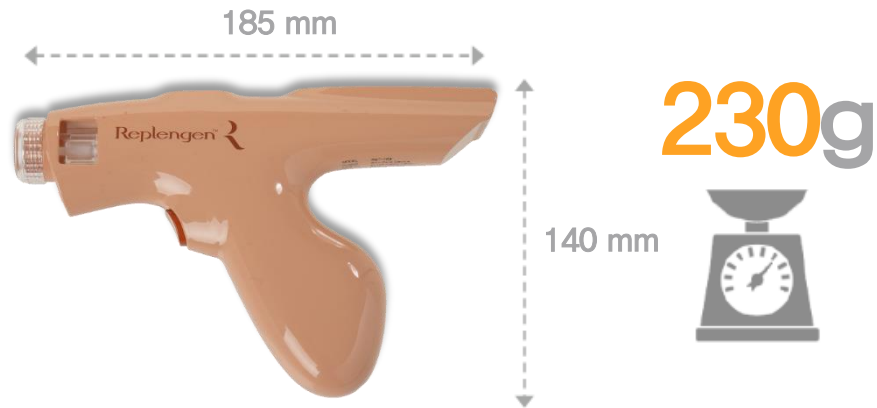
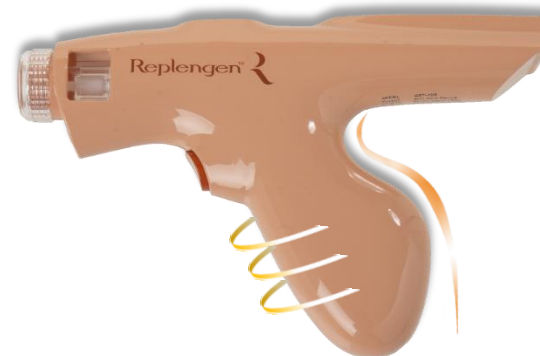
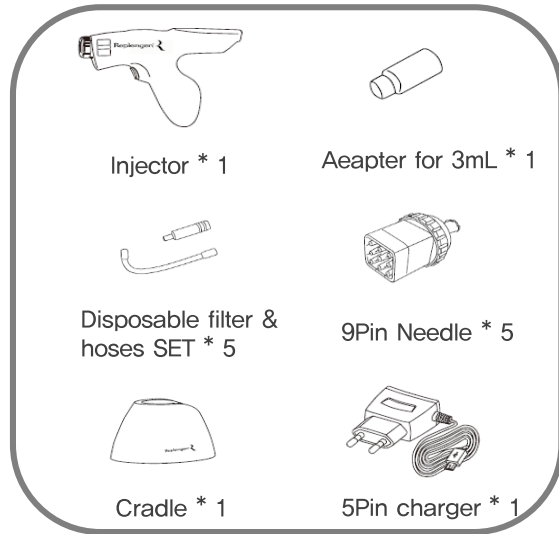
<u>Voltage</u>	▪ 450mmHg/cm2
<u>Charging time</u>	▪ Approx. 60min
<u>Continuous use</u>	▪ Approx. 90min
<u>Applicable syringe</u>	▪ 1mL ~ 3mL

Thank You

Reanzen Co., Ltd.

#503 ~ 505 / 250, Hagui-ro Dongan-gu, Anyang-si, Gyeonggi-do, 14056, South Korea
경기도 안양시 동안구 학의로 250 / 503호 ~ 505호(관양동, 관양두산벤처다임)
Tel : +82. 31. 596. 6310, Fax :+82. 31. 596. 6311, www.reanzen.com